

Designed
with your needs in mind



CrystalClear Plus

Nurses and the Respiratory Care support network in Acute Care

Rüsch CrystalClear Plus Thermosensitive PVC Tracheostomy Tubes have been designed for patient comfort, effective stoma care and additional patient safety features were reviewed at every step - included adding a colour coding feature to allow the tube size to be matched to the correct inner cannula size easily.

Nurses and the Respiratory Care support network have the highest interaction with tracheostomy patients' and their relatives, making it our priority to provide support on the challenges this group faces from a product usability as well as from a daily care routine/treatments perspective.


How were the needs of your profession considered in the CrystalClear Plus product design process?

- All components in the kit are single use to reduce risks of cross infection to patients and carers as no cleaning required.
- Flexible neck flange and a soft neck band reduces movement of the tube against granulating tissue during use, which can help in minimising the risk of bleeding, infection and pain.
- Thin walled cuffed and un-cuffed tracheostomy tubes are available in sizes 3.5mm to 10.0mm and includes all half sizes to support patient comfort and breathing effort.
- Each kit (Size 6 and above) includes 2 flexible inner cannula.

Rüsch CrystalClear Plus Tracheostomy Tube

15 mm Swivel Connector
Designed to reduce torque with ventilator circuits

Cuff
Tubes are available both with or without a cuff

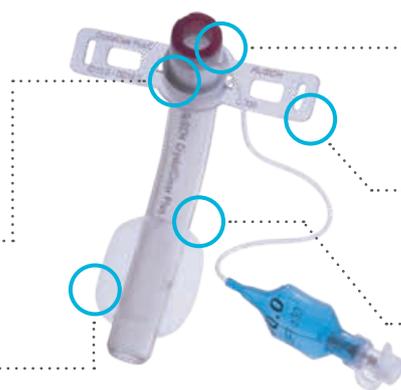

Inner Cannula
Size-specific colour coding to aid size identification

Thin, Flexible Flange
Facilitates effective stoma care

95° Shaft
Promotes a secure fit to help prevent accidental decannulation²

Rüsch CrystalClear Plus Tracheostomy Tube

PRODUCT CODE*	PRODUCT CODE**	SIZE / I.D. (MM)	O.D. (MM)	LENGTH (MM)	ANGLE Ø APPROX.	CUFF Ø (MM)	INNER CANNULA	TUBE	COLOUR CODE	QTY
121300-000035	121303-000035	3.5	5.3	47.0	95°	11.0	-	3.5	Red	5
121300-000040	121303-000040	4.0	5.9	54.0	95°	12.0	-	4.0	White	5
121300-000045	121303-000045	4.5	6.5	54.0	95°	14.0	-	4.5	Grey	5
121300-000050	121303-000050	5.0	7.0	57.0	95°	14.0	-	5.0	Blue	5
121300-000055	121303-000055	5.5	7.5	58.0	95°	15.0	-	5.5	Green	5
121300-000060	121303-000060	6.0	8.7	64.0	95°	19.0	6.0	-	Orange	5
121300-000065	121303-000065	6.5	8.7	64.0	95°	19.0	6.5	-	Maroon	5
121300-000070	121303-000070	7.0	10.0	70.0	95°	21.0	7.0	-	Blue	5
121300-000075	121303-000075	7.5	10.7	75.5	95°	25.0	7.5	-	Yellow	5
121300-000080	121303-000080	8.0	11.3	81.0	95°	28.0	8.0	-	Lilac	5
121300-000085	121303-000085	8.5	11.3	81.0	95°	28.0	8.5	-	Green	5
121300-000090	121303-000090	9.0	12.7	87.0	95°	28.0	9.0	-	Brown	5
121300-000095	121303-000095	9.5	12.7	87.0	95°	28.0	9.0	-	Navy	5
121300-000100	121303-000100	10.0	14.0	95.5	95°	30.0	10.0	-	Red	5

*Without Cuff. **With Cuff.

Kit contains: Tracheostomy Tube (cuffed or un-cuffed, Soft Neck Band, 2 colour coded inner cannula, Patient card, Satchet of water based lubricant gel

Rüsch CrystalClear Accessories

PRODUCT CODE	DESCRIPTION	QTY
121311	Shower cap for Rüsch CrystalClear Plus	4
121312	Cough cap for Rüsch CrystalClear Plus	4
121313	Sealing cap for Rüsch CrystalClear Plus	4

Rüsch CrystalClear Inner Cannula

PRODUCT CODE	SIZE	I.D. (MM)	O.D. (MM)	COLOUR CODE	QTY
121310-000060	6.0	5.0	5.7	Orange	10
121310-000065	6.5	5.0	5.7	Maroon	10
121310-000070	7.0	6.0	6.7	Blue	10
121310-000075	7.5	6.5	7.3	Yellow	10
121310-000080	8.0	7.0	7.7	Lilac	10
121310-000085	8.5	7.0	7.7	Green	10
121310-000090	9.0	8.0	8.7	Brown	10
121310-000095	9.5	8.0	8.7	Navy	10
121310-000100	10.0	8.8	9.7	Red	10

www.teleflextracheostomy.com



References

1. Wilkinson KA, Martin IC, Freeth H, Kelly K, Mason M. NCEPOD: On the right Trach? 2014. Available at: <http://bit.ly/ontherighttrach>. Accessed August 22, 2017.
2. Mallett J, Albarran JW, Richardson A, eds. Critical Care Manual of Clinical Procedures and Competencies. Oxford, UK: Wiley Blackwell; 2013.

Please ask your Teleflex Sales Representative if you wish to know more or visit our Rüsch Tracheostomy micro site where you will find out more about this range and our other Rüsch exciting products. Contents of unopened, undamaged package are clean, ready to use and for single-use. Refer to package insert for current warnings, indications, contraindications, precautions, and Instructions for Use.

Teleflex, the Teleflex logo, CrystalClear and Rüsch are trademarks or registered trademarks of Teleflex Incorporated or its affiliates, in the U.S. and/or other countries. Information in this document is not a substitute for the product Instructions for Use. The products in this catalogue may not be available in all countries. Please contact your local representative. Subject to technical changes without further notice. Revised 02 2023

© 2023 Teleflex Incorporated. All rights reserved. MCI-101150 · REV 0 · 02 23 PDF

Teleflex

Distributed by:
Teleflex Headquarters International, Ireland · Teleflex Medical Europe Ltd. · IDA Business and Technology Park
Dublin Road · Athlone · Co Westmeath · Tel. +353 (0)9 06 46 08 00 · Fax +353 (0)14 37 07 73 · orders.intl@teleflex.com
United Kingdom Tel. +44 (0)14 94 53 27 61 · info.uk@teleflex.com
South Africa Tel. +27 (0)11 807 4887 · assist.africa@teleflex.com