

Designed
with your needs in mind



TracFlex Plus


Clinicians - Acute Care units, ENT, OR

Teleflex have continued innovating and improving our designs over the past 35 years, to now offer our newly upgraded range of TracFlex Plus Adult cuffed and un-cuffed tracheostomy tube sets.

Clinical indications for tracheostomy tube use is diverse and patient needs often change, we design our products to achieve improved patient outcome. Our belief is that patient's comfort and acceptance of their tube not only helps in their recovery but also supports your teams in the daily care of patients, supporting relatives and reducing the risk of complications, making the whole experience easier for all.

How were the needs of your profession considered in the TracFlex Plus product design process?

- The tube is easily controlled during insertion - has a smooth conical tube tip and a guide wire channel in the obturator.
- Flexible/non kink tubing assists in patient comfort and tube acceptance for both usual and unusual anatomies.¹
- Pliable neck flange reduces movement of the tube against granulating tissue during use, which can help in minimising the risk of bleeding, infection and pain.²
- The adjustable flange moves freely when unlocked for smooth positioning and an audible click confirms when securely locked.

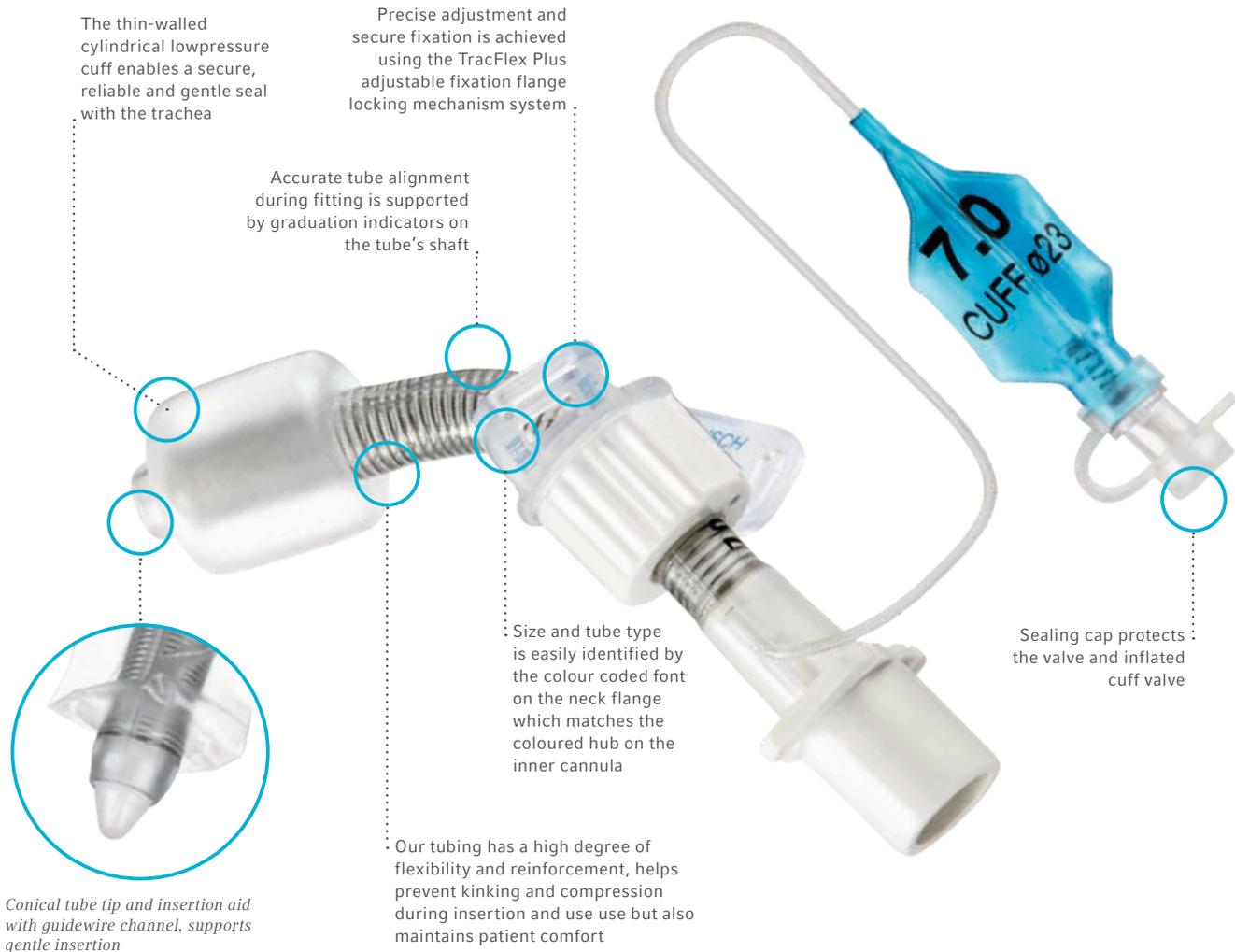

TracFlex Plus Set (with Cuff, Ref. 121903; without Cuff, Ref. 121900)

ORDER SIZE / I. D.	O.D. IN MM	LENGTH (CL) RANGE IN MM	DIM. A IN MM	DIM. B IN MM	DIM. C IN MM	ANGLE θ (APPROX.)	CUFF Ø IN MM **
● 7.0	10.1	60.3 - 66.5	37.5	6.2	22.8	100°	23
● 8.0	11.1	66.9 - 88.5	41.5	21.6	25.4	100°	27
● 9.0	12.1	70.4 - 117.3	44.5	46.4	26.4	100°	30
● 10.0	13.1	78.5 - 121.3	44.5	42.8	34.0	100°	32
● 11.0	14.2	80.6 - 125.3	47.5	44.7	33.1	100°	34

** Only cuffed versions

TracFlex Plus Inner Cannulas (Ref. 121910)

ORDER SIZE	I.D. IN MM	O.D. IN MM	QTY
● 6.0	5.0	5.7	10
● 7.0	5.5	6.2	10
● 8.0	6.5	7.2	10
● 9.0	7.8	8.3	10
● 10.0	8.8	9.4	10
● 11.0	9.8	10.4	10

Our Flexible Inner Cannula System: Effective Secretion Management Made Easy

TracFlex Plus inner cannula is designed to protect the inner lumen of the tracheostomy tube against the build up of secretions. Single use to support good hygiene, reducing the risk of infection to patients and carers.

The TracFlex Plus thin walled, flexible, inner cannula enables a close fit to the inner diameter of the tracheostomy tube, with minimal torque on the tube during insertion and removal of the inner cannula.

TracFlex Plus Accessories

REF.	DESCRIPTION	QTY
121911	Shower cap for TracFlex Plus tubes, acting as a splash guard to prevent the entry or inhalation of water droplets during bathing or showering	4
121912	Cough cap for TracFlex Plus tubes, to channel expelled secretions downwards	4
121913	Sealing cap for TracFlex Plus tubes, for weaning therapy	4

Please ask your Teleflex Sales Representative if you wish to know more or visit our Rusch Tracheostomy micro site where you will find out more about this range and our other Rusch exciting products.

Teleflex, the Teleflex logo, Rusch, and TracFlex Plus are trademarks or registered trademarks of Teleflex Incorporated or its affiliates, in the U.S. and/or other countries. Information in this document is not a substitute for the product Instructions for Use. The products in this catalogue may not be available in all countries. Please contact your local representative. Revised 12 22. Subject to technical changes without further notice. © 2023 Teleflex Incorporated. All rights reserved. MCI-101148 · REV 0 · 02 23 PDF